


FAQ

OB Application

**Why is the Auto Annotation Button
not Visible on the Touchscreen**

[FAQ – Auto Annotation in OB]

- Question
I don't see the **Auto Annotation** button on the touchscreen. What should I check?
- Answer
 - ▶ The button appears in **OB Application** in **Freeze** mode or while reviewing saved images.
 - Go to **Utility > Touch Customization** and check if the **Auto Annotation** button is assigned.
 - Go to **Setup > Assist > OB > General** and make sure **ViewAssist** is set to **On**.

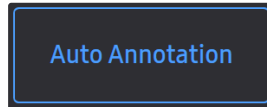

- If **ViewAssist™** is already **On** but the button is still not visible, check the **Auto Annotation** tab and ensure the checkboxes for the desired views are selected.
- You can also press the **Edit** button to modify the **View Title** or **Segmentation Label** for each view.